



Comprehensive  
and Intuitively Designed

# Gallery™ Laminoplasty Spine System



# Gallery Laminoplasty®

## Spine System

The Gallery Laminoplasty Spine System is a comprehensive and intuitive system designed to prevent spinal cord impingement in the lower cervical and upper thoracic spine.



### Flexibility

- Cobra head design
- Variety of plate styles and sizes
- Self-tapping and self-drilling screws




### Simplicity

- Design incorporates quick, secure connection between the screws and inserter
- Spin Tight Technology makes screw insertion easy

## Smart Plate Design

The Gallery Laminoplasty Spine System offers two plate styles. Both cobra head plate designs are low profile and offer multiple points of fixation for ideal screw placement.

- Design incorporates a quick, secure connection between the screws and inserter
- Spin Tight Technology makes screw insertion easy

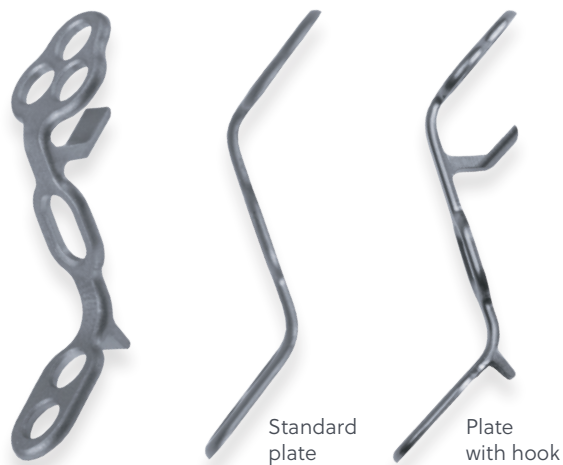


## Variety of Plate Sizes and Styles

The Gallery Laminoplasty Spine System plates are available in hook or standard configurations, in a variety of lengths from 8 mm to 16 mm (2-mm increments).




## Self-Tapping and Self-Drilling Screws

Gallery Laminoplasty Spine System screws are available in self-drilling and self-tapping options in multiple lengths from 5 mm–11 mm (2-mm increments). Screws are available in 2.4 mm and 2.8 mm diameters.

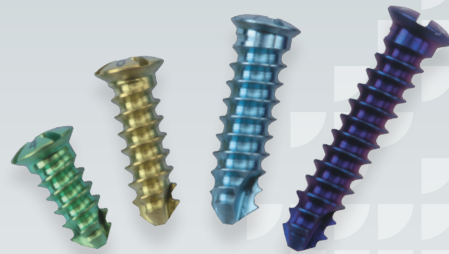

Self-Tapping Screws

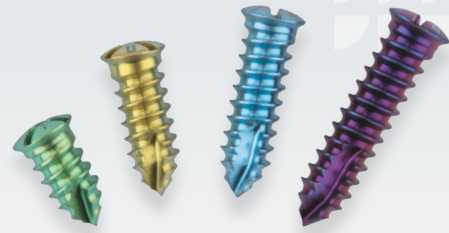

Self-Drilling Screws



## Ergonomic Inserter

### Quick, Secure Connection

The design behind the Gallery Laminoplasty Spine System inserter ensures a quick, secure interface between the screw and inserter.

- Quick, secure connection
- Wheel mechanism
- Ergonomic handle fits in the palm of your hand

### Spin Tight Technology

Spin Tight Technology offers a wheel mechanism on the handle, the location of which allows the surgeon to simultaneously place pressure on the inserter with the palm of the hand and drive the screw with minimal effort.

For more information, visit [ZimVie.com](https://www.ZimVie.com)

**ZimVie**  
10225 Westmoor Drive  
Westminster, CO 80021



All content herein is protected by copyright, trademarks and other intellectual property rights owned by or licensed to Zimmer Biomet Spine, Inc. (d/b/a ZimVie Inc.) or one of its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet Spine, Inc. (d/b/a ZimVie Inc.). This material is intended for health care professionals, the ZimVie Spine sales force, and authorized representatives. Distribution to any other recipient is prohibited.

For product information, including indications, contraindications, warnings, precautions, potential adverse effects, and patient counseling information, see the package insert and [www.ZimVie.com](https://www.ZimVie.com). ZV0828 REV A 03/23 ©2023 ZimVie Inc. All rights reserved.